



### Polypropylene barrel pumps BÜRKLE

Suitable for acids, basis and detergents transfers. With rigid discharge tube. Container adaptor not included.

Submersible length (cm)	Total length (cm)	Pumping cap. (ml/line)	Pack (u.)	Art. No.
100	109	300	1	5600-01001
50	59,5	200	1	05600-0501

#### Required adapter to choose



PVC adapter



2" internal thread adapter




R2" internal thread adapter



Trisure internal thread adapter

Description	Pack (u.)	Art. No.
PVC barrel plug for openings with a diameter from 40-70mm	1	0005600-15
Mauser 2", internal coarse thread	1	0005600-17
R2", internal steel barrel thread	1	0005600-13
Trisure, internal coarse thread	1	0005600-18

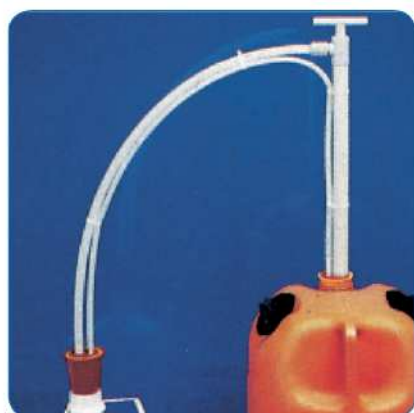

### Polypropylene blocked pump to prevent vapour escapes BÜRKLE

Polypropylene, with a PVC hose. Evacuation gases created during the pumping process are returned to the original container through a gas pipe so that harmful or smelly gases do not escape outside. With a universal adaptor for openings from 30 to 70mm. Suitable for containers from 40 to 70cm in height. Hose length 120cm approximately. Recommended by BG Chemie Germany.

Description	Pack (u.)	Art. No.
PP blocked pump	1	004-0E0131



### Pump-it® container pump BÜRKLE

The manual pumping of the cylindrical plunger allows air to enter and liquid extraction by positive pressure. Materials in contact with liquid ; PP, PE and Viton®. For liquids with viscosity up to 1500cps at 20°C. 64mm diameter thread (BSI). Suitable for weak acid, base and hydrocarbon transfers. Pumping capacity: 8L/min. Immersion tube of 100cm. Adaptors available for different containers and standard canisters.

Description	Pack (u.)	Art. No.
Pump-it® pump	1	004-016011

→ Continued on the next page

For further information, please contact [scharlab@scharlab.com](mailto:scharlab@scharlab.com)

**Adapters for outer thread of container**


Description	Pack (u.)	Art. No.
2" thread, thick Mauser L-ring type (green)	1	004-016041
2" thread, Van Leer Trisure type (brown)	1	004-016061
2" thin thread, for metal drums (blue)	1	004-016031

**Adapters for inner thread of container**


Description	Pack (u.)	Art. No.
Thread DIN60, Mauser Faset (black). Suitable for 343-025NMA (left)	1	004-016051
Thread DIN51 (black) (right)	1	004-0H0711
Thread DIN61 (yellow/blue)	1	0005201-21
Thread DIN71 (brown/blue)	1	0005201-22
Adapter set	5	004-016021


**Stainless steel pumps for flammable liquids  
BÜRKLE**

Pumps for steel containers up to 200L.

Description	Submersible length (cm)	Total length (cm)	Pumping cap. (ml/line)	Pack (u.)	Art. No.
Steel pump with rigid spout (left)	91	103	560	1	0005601-80
Steel pump with flexible 1,2m PTFE/steel tube and NW8 ball valve with a 12,8mm outlet (right)	91	103	560	1	0005601-081
Required washer, for use with steel pump with flexible PTFE/steel tube	-	-	-	1	05601-0100

 Visit our new web [www.scharlab.com](http://www.scharlab.com)